



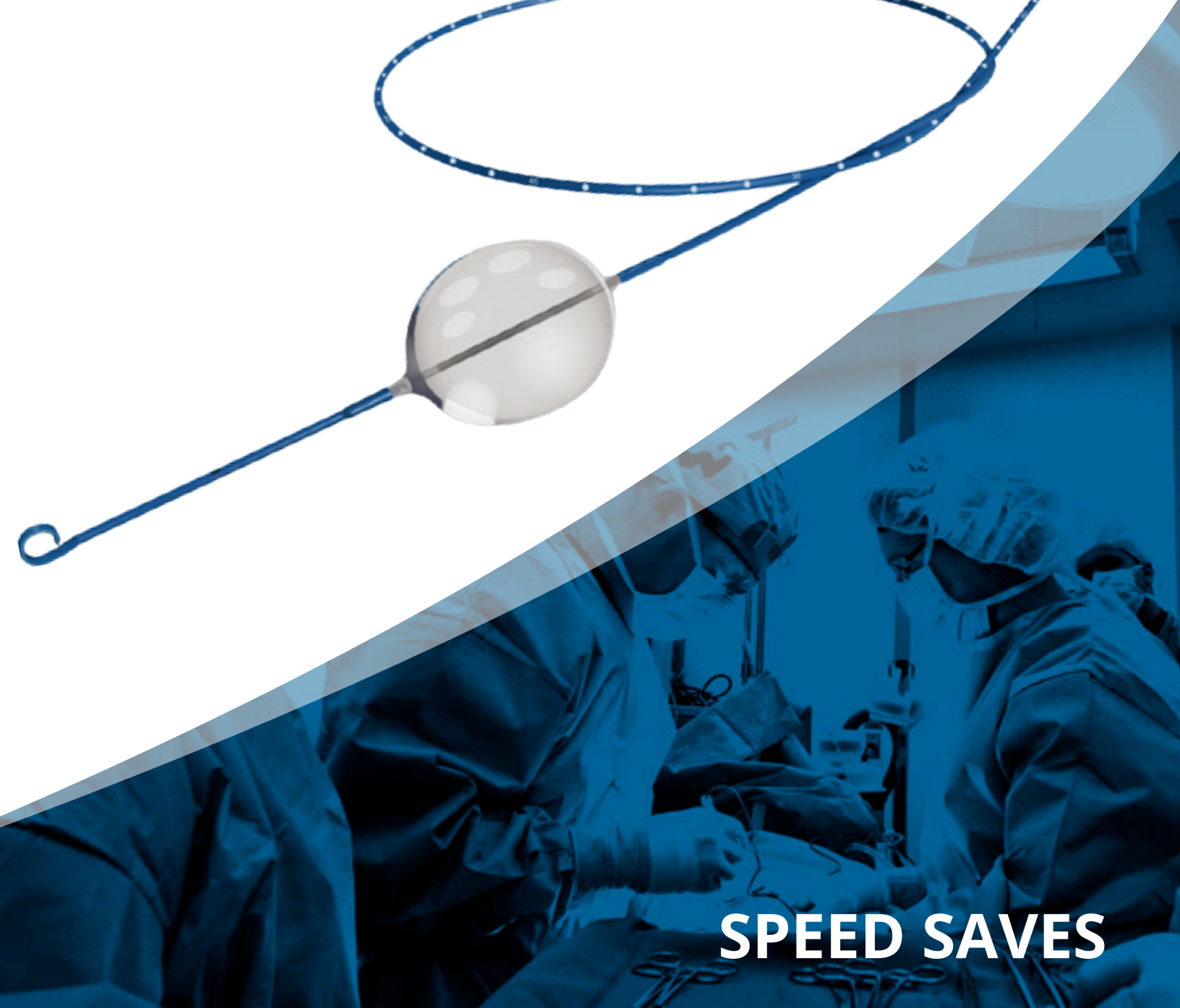

PRYTIME
MEDICAL™

The REBOA Company™

ER-REBOA™ Catheter

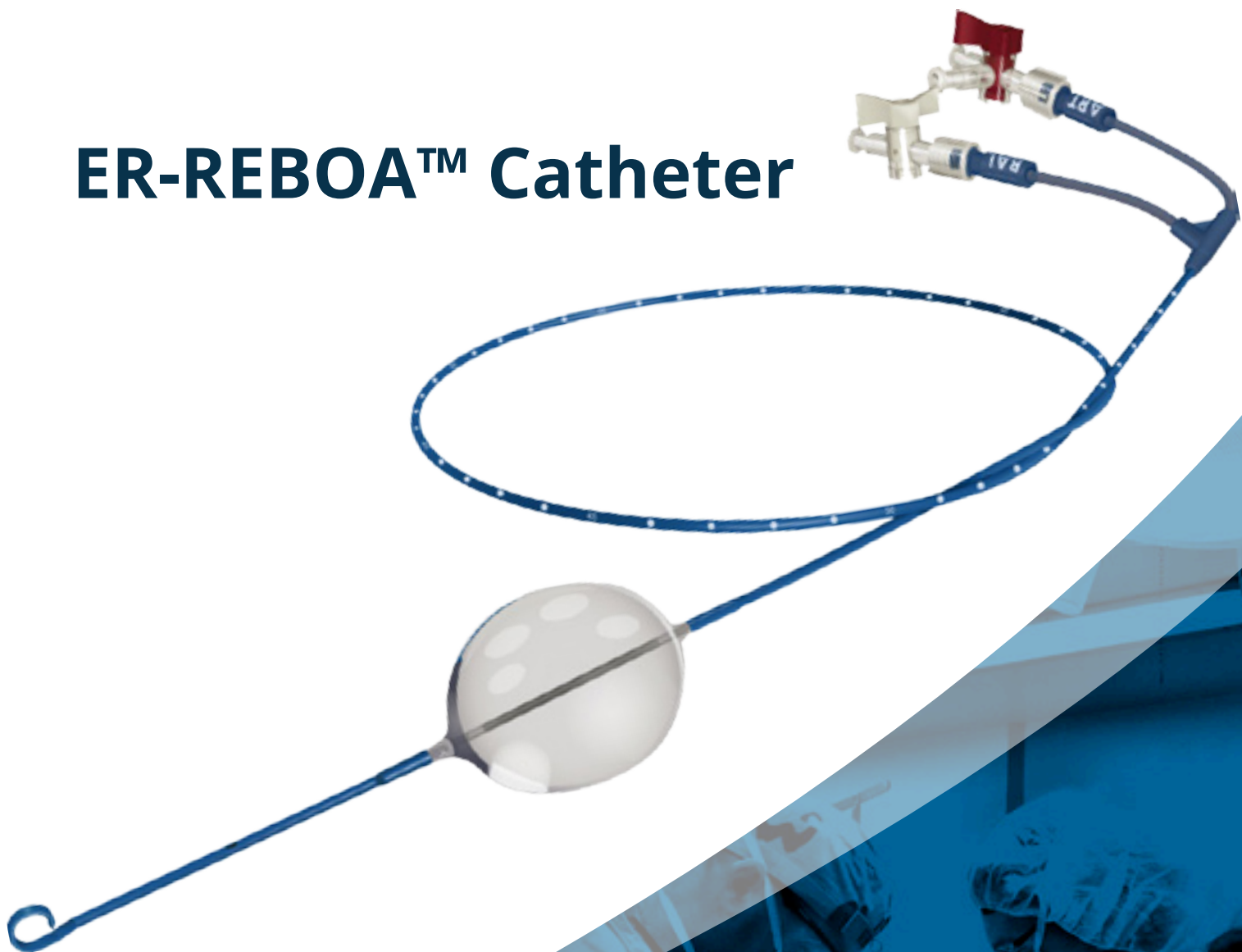


SPEED SAVES

REBOA

Resuscitative Endovascular Balloon Occlusion of the Aorta

A minimally invasive technique using a balloon catheter to temporarily occlude large vessels and control hemorrhage for patients in traumatic arrest and hemorrhagic shock states.

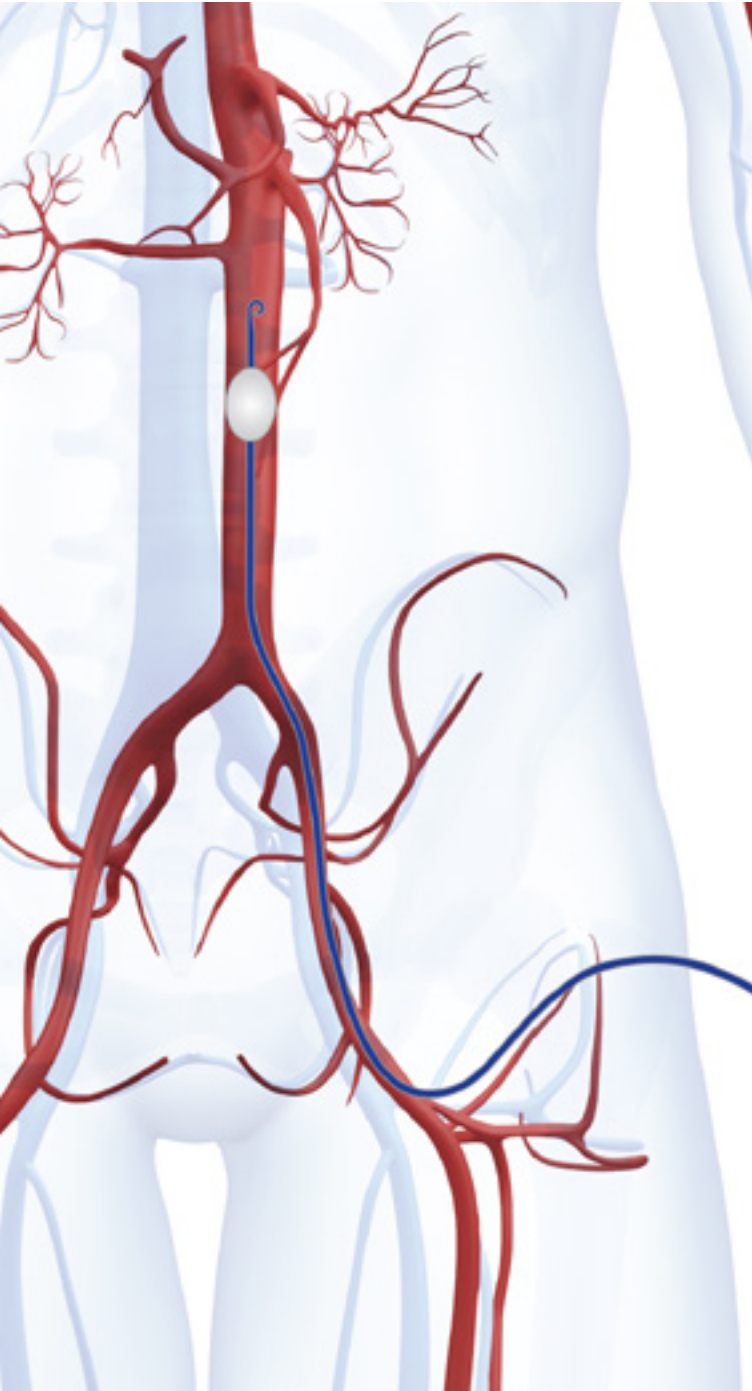

Balloon Landing Zones

Aortic Zone 1

Extends from the origin of the left subclavian artery to the celiac artery (approximate vessel diameter 20mm – young adult male)¹

Aortic Zone 3


Extends from the lowest renal artery to the aortic bifurcation (approximate vessel diameter 15mm - young adult male)¹

ER-REBOA™ Catheter

➔ Designed for use in the emergency and critical care environments

➔ Optimized to support rapid and immediate hemorrhage control

- 7Fr sheath compatibility – Reduces the need for surgical repair at the access site.²
- Guidewire free design – Easy to use; facilitates rapid placement in emergent situations.
- Atraumatic P-tip® – Prevents catheter migration into branch vessels, and aids in correct positioning to desired location.
- Compliant balloon – Material is designed to support inflation while minimizing injury to the vessel.
- Length marks – Allow physicians to advance the catheter to the desired position.

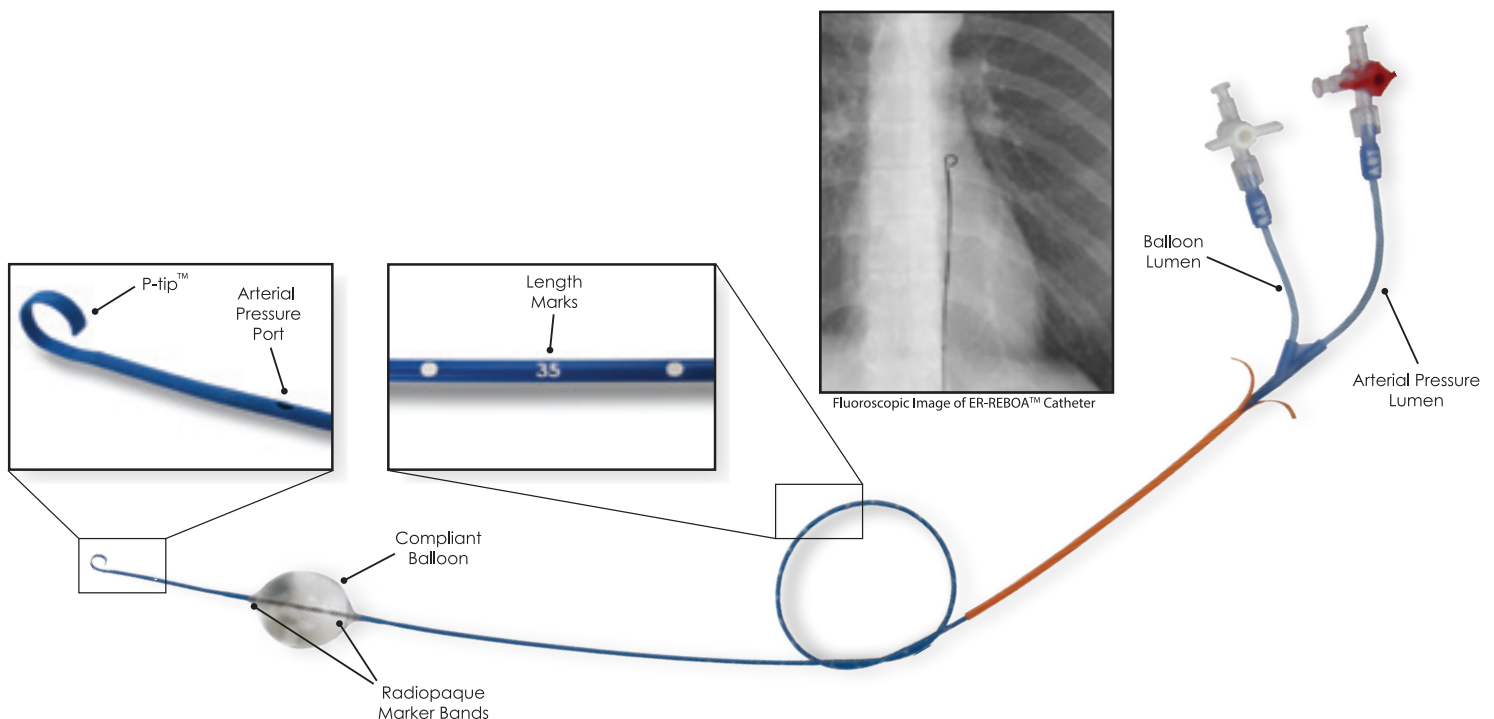



No one should bleed to death.

Ordering Information

ER-REBOA™ Catheter Order Number #ER7232A

Convenience Kit Order Number # KT1835C

sales@prytimemedical.com

Phone: 210-340-0116

www.prytimemedical.com

